

HEADQUARTERS

Apex Medical Corp.
No. 9, Min Sheng St., Tu-Cheng,
New Taipei City, 23679, Taiwan
T +886 2 2268 5568
F +886 2 2268 6525
marketing@apexmedicalcorp.com

PAN-ASIA

TAIWAN
Apex Medical Corp.
No. 9, Min Sheng St., Tu-Cheng,
New Taipei City, 23679, Taiwan
T +886 2 2268 5600 /0800 054 133
F +886 2 2268 6526
customer@apexmedicalcorp.com

STURDY INDUSTRIAL
168, Sec.1, Zhongxing Road, Wugu District,
New Taipei City 24872, Taiwan
T +886 2 2981 5245
F +886 2 2984 8208
foreign@sturdy.com.tw
www.sturdy.com.tw

CHINA
Apex Medical (Kunshan) Corp.
No. 1368, ZiZhu Road, Kunshan,
Jiang Su Province, China
T +86 512 8617 1660 Ext 1150
F +86 512 8617 8663
customer.cn@apexmedicalcorp.com

INDIA
Apex Medicalcorp India Pvt. Ltd.
303, Third Floor Block 4A DLF
Corporate Park Phase-III,
Gurgaon, Haryana
T +91 11 6654 4138
F +91 11 6654 4052
enquiries.india@apexmedicalcorp.com

THAILAND
Apex Medical Thailand
65/100 Chamnanphenjati Business Center Bldg.
11/F Rama 9 Road, Huaykwang
Bangkok 10320
T +662 948 4848
F +662 948 4170
sales.thai@apexmedicalcorp.com

MIDDLE EAST • AFRICA
29, Misr Liltameer Building,
Stree#10, AREA#8, Sheraton Building
Heliopolis Cairo, Egypt.
T +202 2 2671303
F +202 2 2678649
raafat.wahby@apexmedicalcorp.com

PAN-EUROPE

SPAIN
Apex Medical S.L.
Elcano 9, 6ª planta 48008
Bilbao. Vizcaya. Spain
T +34 94 4706408
F +34 94 4706409
info@apexmedial.es


UK
Apex Medical Limited
Unit 33, Great Western Business Park
McKenzie Way, Worcester ,WR4 9GN
United Kingdom
T +44 (0) 190 5774 695
sales@apexmedicalcorp.co.uk

FRANCE
Apex Medical France
4 Boulevard de la Chanterie
49124 Saint-Barthelemy d'Anjou
France
T +33 2 72 79 27 13
F +33 2 85 73 07 67
info.france@apexmedicalcorp.com

GERMANY
SLK Vertriebsgesellschaft GmbH
Lindenhorster Straße 38-40 44147
Dortmund
T +49 (0)231-925360-175
F +49 (0)231-925360-29
info@slk-gmbh.de
www.slk-gmbh.de

AMERICAS

USA
Apex Medical USA
927 Mariner Street,
Brea, CA 92821, USA
T +1 714 671 3818
F +1 714 494 8173
sales.usa@apexmedicalcorp.com

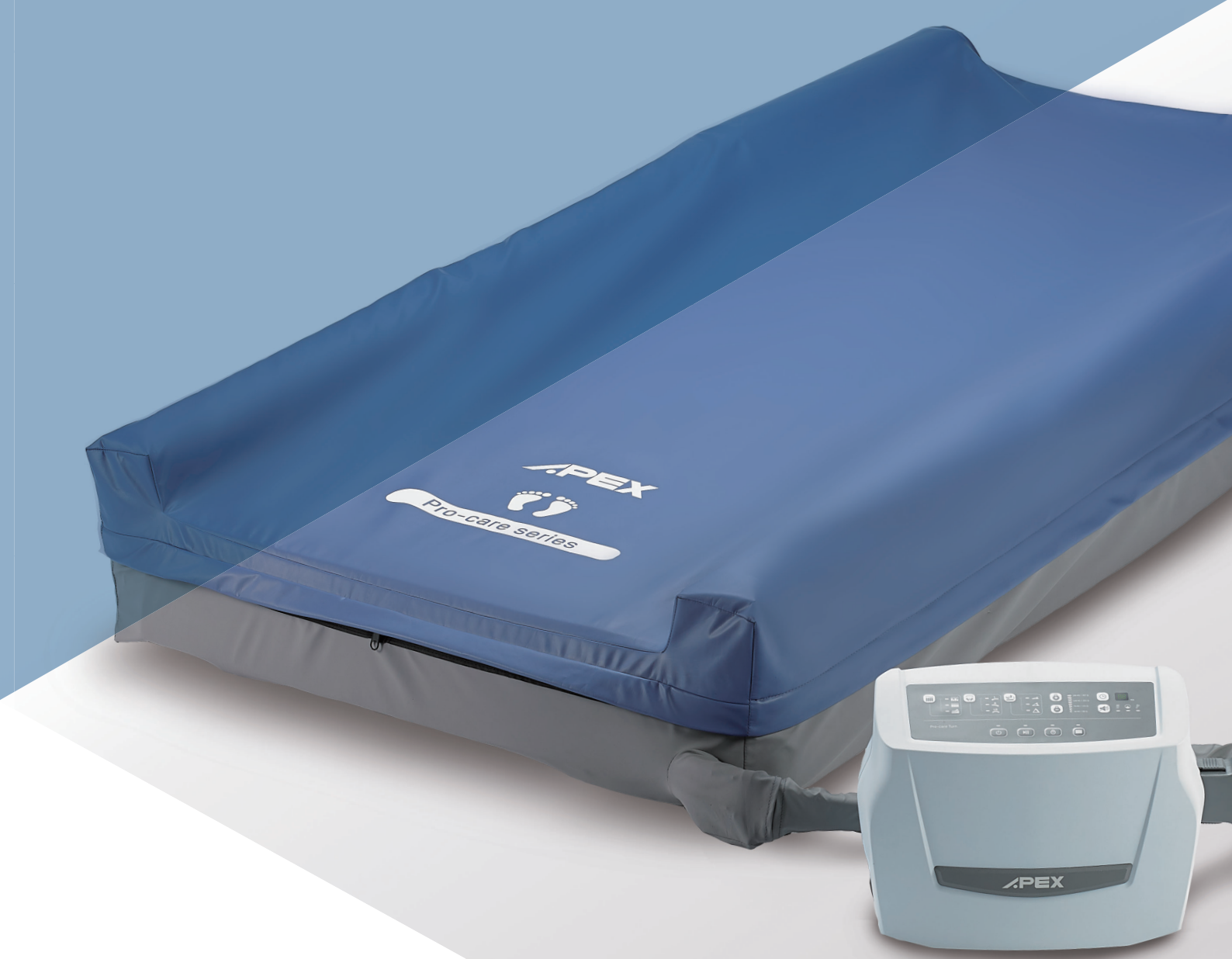

Pro-care Turn

Bilateral Turning Air Mattress

- Suitable for patients in high to very high risk of pressure injuries
- 30° of bilateral rotation therapy to enhance pressure relief
- 3-in-1 therapy available for individual healthcare plans with effortless patient repositioning
- Designed for use in the acute sector



Care for a Healthy Life






Pro-care Turn

Bilateral Turning Air Mattress

“A Turning Point for Better Patient Care”

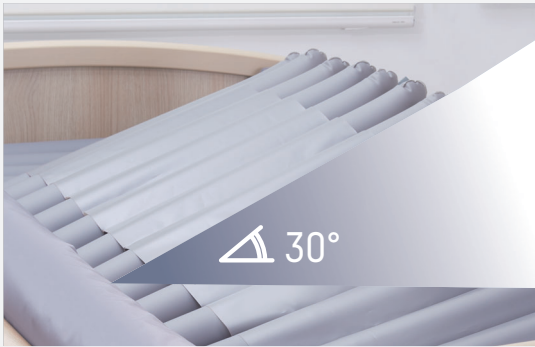

Enhancing Pressure Relief by 30° of Bilateral Rotation

Bilateral rotation can significantly reduce interface pressure on patients’ body whilst increase capillary circulation. ^{*1*2}
30 degrees of turning not only complies with clinical guidelines but also provides optimal comfort. ^{*3}

References:
^{*1} Peterson MJ et al., Effects of turning on skin-bed interface pressures in healthy adults. Journal of advanced nursing. 2010
^{*2} Peterson MJ et al., Patient repositioning and pressure ulcer risk—Monitoring interface pressures of at-risk patients, JRRD 2013
^{*3} Prevention and Treatment of Pressure Ulcer: Clinical Practice Guideline (NPUAP/EPUAP/PPPIA 2014)

4-way Stretch Cover

- Antimicrobial
- Low friction and shear
- Vapour permeable, water resistant
- High-flame retardant

Antimicrobial

4-Way Stretch

Vapour Permeable

Flame Retardant

Effective Pressure Relief by Bilateral Rotation
The pressure mapping below shows a clinical proven body pressure relief by 30° of turning

Position	Left Turn	Supine	Right Turn
<32mmHg (30°)	97%-99%	95%	97%-99%
Pressure Mapping (Weight 75 Kg)			

Seat Inflation

Offering extra support in sacral area during sitting position

Cable Management

Reducing tripping hazards for a safer healthcare environment

Comfort Stability

The Side Bolsters

Thorough patient protection by reducing the chance of being trapped or compressed by the bed rails. Compatible with profiling beds

The Rectangular Cells

Offering comfortable pressure relief by evenly redistribute body pressure

Heel Relief Function

Single air-cell deflation can float the heels with **zero pressure**

Visual and Audible Alarm

- Low Pressure
- Power Failure
- Tech Support

Transport Mode

Non-stop static pressure relief for **more than 6 hours** during power failure or transport



+ +

Effortless, 3-in-1 Therapy Solutions

- Alternating, turning, continuous low pressure are available to tailor to individual healthcare plans

Customised for Optimal Performance

- Selectable turning angles and pressure levels for individual care plan
- Allowing caregiver to facilitate repositioning easily

Rotational Pause Function

- In any scenario possible, it allows caregiver to stop a turn right away

Skin Safety at All Times

- Alternating mode automatically initiated to ensure non-stop pressure relief during rotational therapy when body shifted to one side

Specification

Model	Pro-care Turn	5" Replacement
Pump	Dimensions	34.1 x 16.5 x 26.0 cm
	Weight	4.5 kg / 9.9 lbs
	Case Material	Flame retardant ABS
	Supply Voltage	230V / 50Hz
	Operating Cycle Time	10, 15, 20, 25 min.
	Dwell Time	10, 20, 30, 60 min.
Mattress	Dimensions	200 x 90 x 13 cm, side bolster 25 cm
	Cell Height	13 cm cells with turning cells on the bottom from 4 th to 12 th
	Weight	11.4 Kg / 25 lbs
	Cover material	4-way stretch
	Cell Material	TPU
	Maximum Patient Weight	180 kg / 397 lbs
	Flame Retardant Standards	EN 597-1 & EN 597-2

Specifications are subject to change for improvement without notice